

Fig.1.

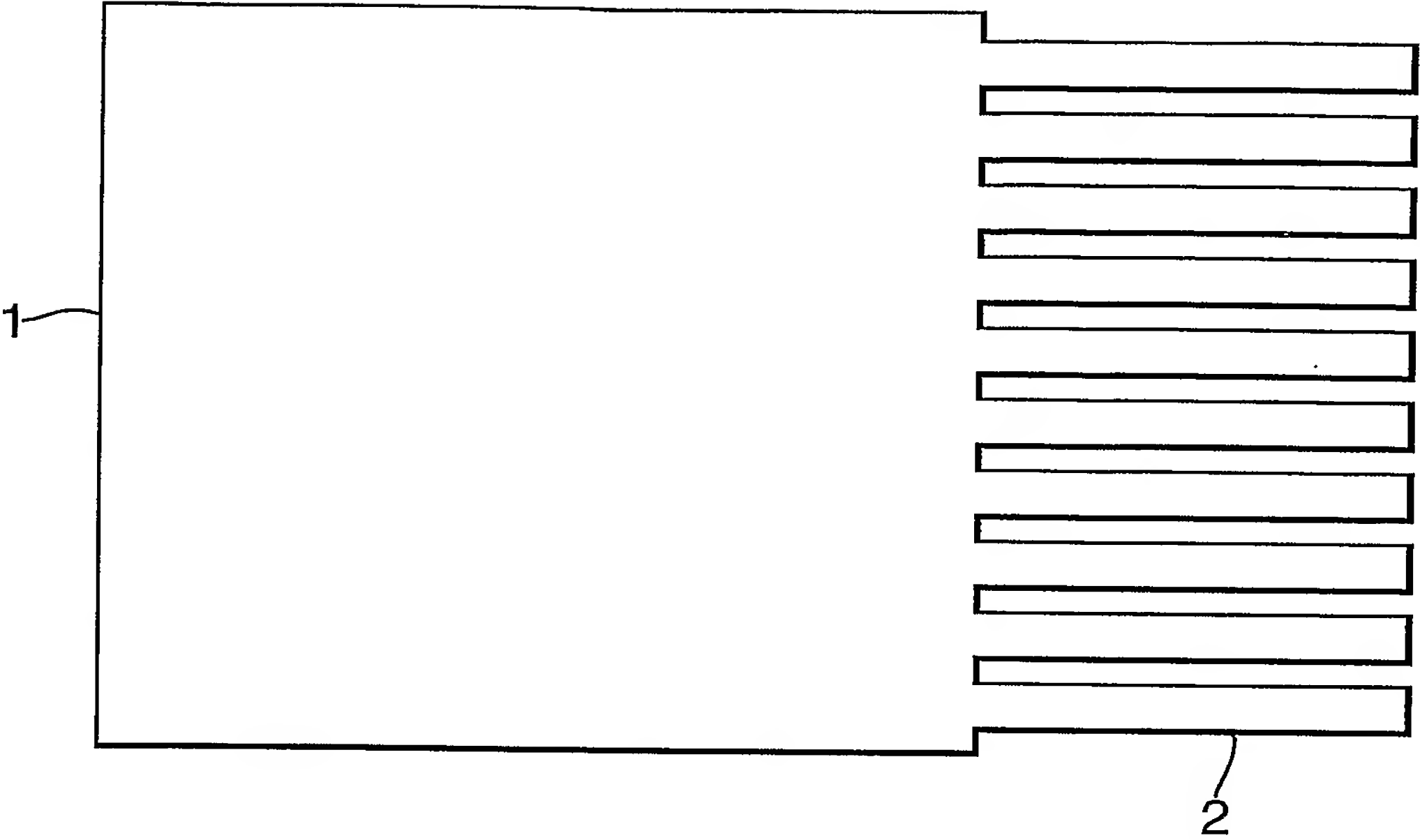


Fig.2.


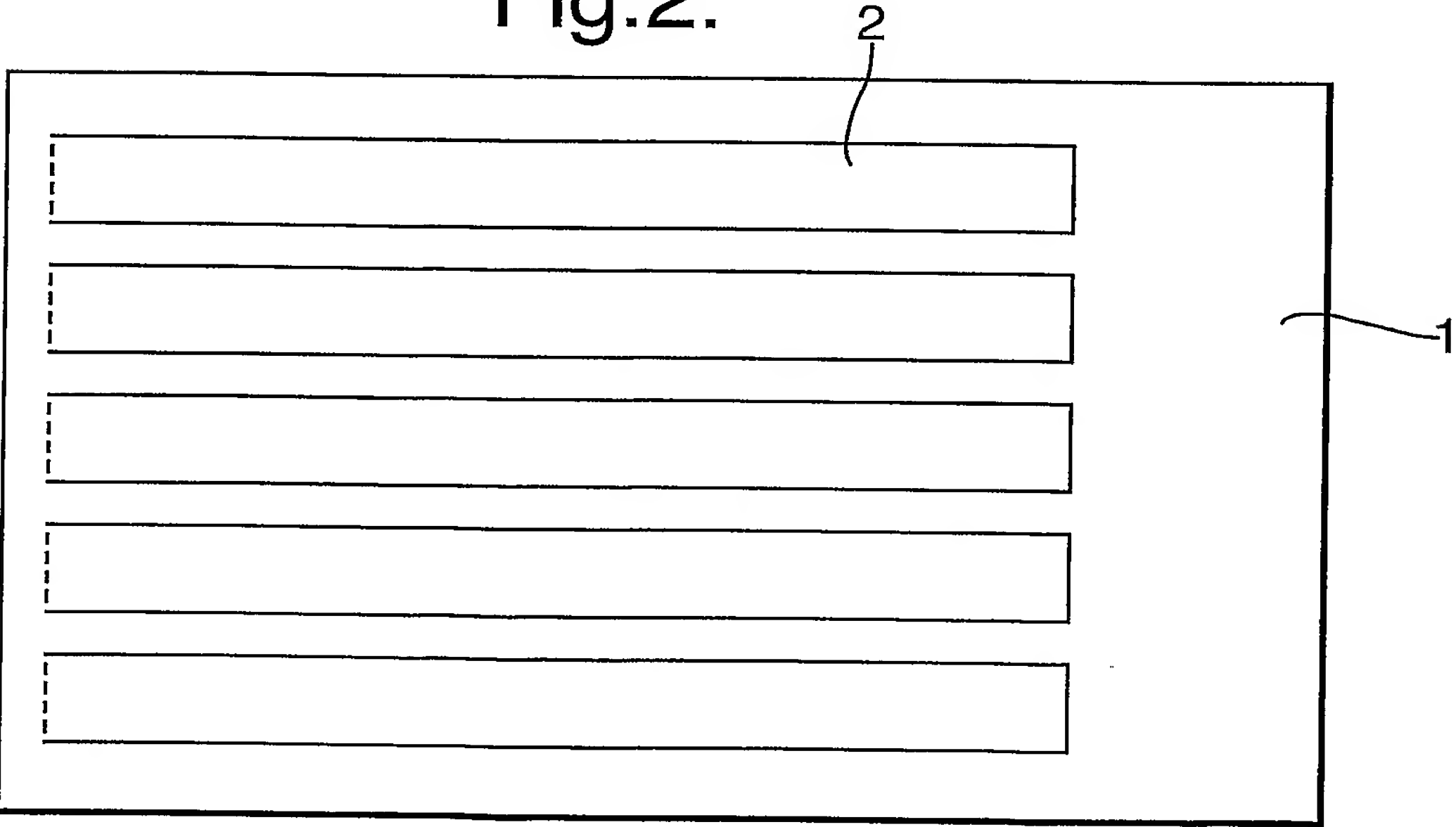


Fig.3.

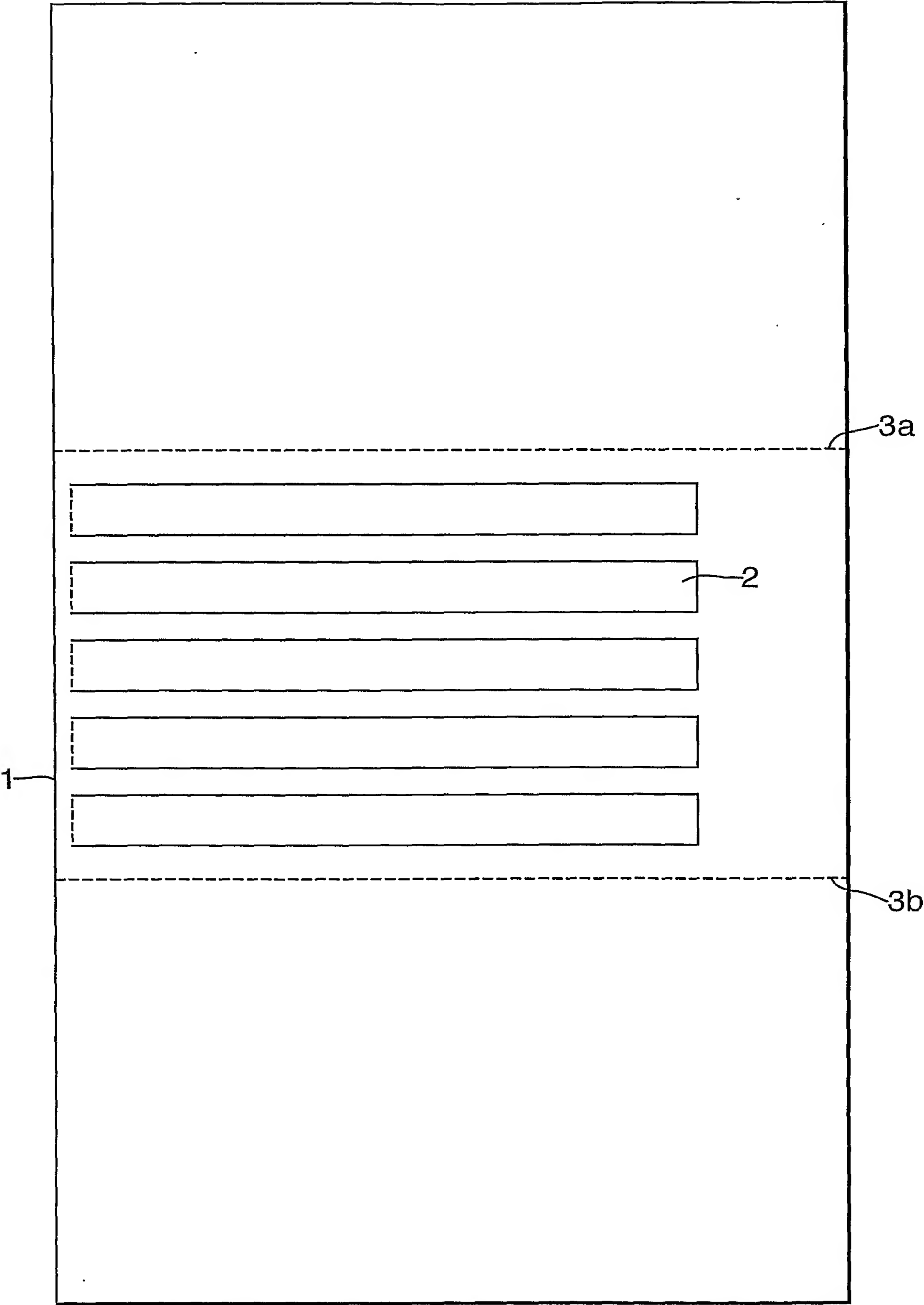

3/4

Fig.4.

WHAT THE DRINK DETECTOR DOES AND ITS LIMITATIONS
 The Drink Detector is a fast, accurate, standard detector of benzodiazepines, ketamine, and GHB. It is a simple, easy-to-use, and reliable device that can be used in a variety of settings, including law enforcement, medical, and forensic laboratories.

Limitations: The Drink Detector is not a definitive test for benzodiazepines, ketamine, and GHB. It is a screening test that can indicate the presence of these substances in a drink. A positive result indicates that the drink contains one or more of these substances, while a negative result indicates that it does not. The Drink Detector is not designed to detect other substances, such as alcohol, opiates, or stimulants.

Instructions for Use:

Benzo Test
 Tear the strip from the card and insert the white end of the strip into the beverage as far as the maximum level indicator (solid blue line) until the purple dye shows on the middle band of the test. Remove from the drink and lay flat until dye has cleared and read results at 5 minutes. Two purple lines shows negative, one line shows positive.

Positive ☐ Negative ☒ Invalid ☐


Ketamine Test
 Dab a spot of the beverage onto the paper and observe the result. A colour change to orange is considered a positive.

Positive ☐ Negative ☒

GHB Test
 Dab a spot of the beverage onto the paper and observe the result. A colour change to blue is considered a positive.

Positive ☐ Negative ☒

Fig.5.

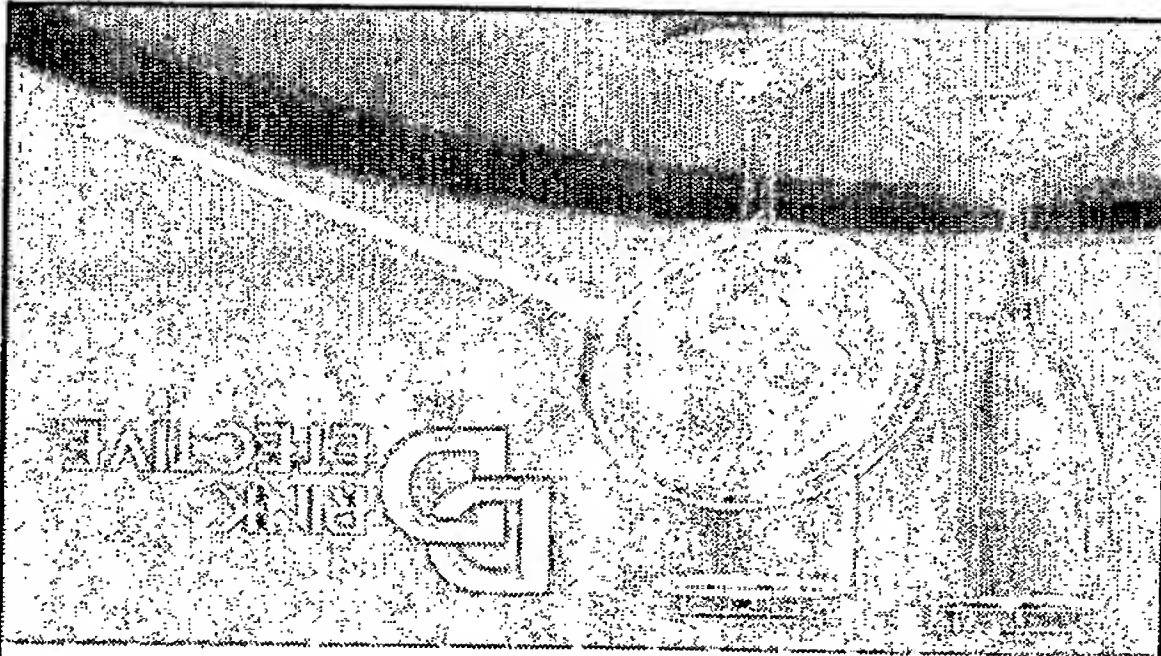

IF YOU HAVE A POSITIVE TEST RESULT

- DO NOT DRINK THE BEVERAGE. TRUST NO-ONE
- Give the drink to someone responsible. Consider phoning the police.
- A drug can take effect in as little as 20 minutes. If you drank from a spiked drink, ask a trusted friend to seek medical advice.
- If you are not with a trusted friend, tell the bar Manager and seek medical advice.
- IF YOU THINK YOU HAVE A SPIKED DRINK, YOUR ONLY CONCERN SHOULD BE YOUR PERSONAL SAFETY.

Sales help to support The Roofie Foundation's work with victims of drug-assisted assault. Visit www.roofie.com www.drinkdetective.com or call our helpline: +44(0)1332 362832

© B.I.L. 2004

CAUTION

Instant Blue colour change

Put 1 drop of the drink onto the centre of the test pad.

HOW TO TEST

HOW TO TEST **K TEST**

Put 1 drop of the drink onto the centre of the test pad.

RESULTS

Red/Orange colour change
May have white centre

CAUTION

Test here

1 drop

HOW TO TEST **B TEST**

Using the dropper, put 3 drops onto the test area as shown below (for best results leave 5 seconds between drops). After 2 minutes read against results shown below.

RESULTS

2 lines - OK

1 line - CAUTION

Test here

3 drops

WHAT THE DRINK DETECTIVE(TM) DOES AND ITS LIMITATIONS

The Drink Detective(TM) is designed to test for the presence of three substances: benzodiazepines (which include Rohypnol and Valium), illicit GHB (gamma-hydroxybutyrate) and ketamine. The Drink Detective(TM) does not test for any other drugs or substances. It is designed to test for minimal standard dosages of these substances and quantities below these levels may not be detected. Due to the sensitivity of the test, it is possible that the presence of certain substances could produce a false positive. In particular, orange and tomato juices and certain milk products might produce a false positive for ketamine. Milk products and soda waters containing bicarbonate of soda might produce a false positive for GHB. Test your favourite drink and be alerted to different reactions. Using this test does not remove your responsibility to drink safely. Follow instructions exactly. Do not alter the test. Do not place in a drink or eat the test. Keep the test away from children, animals, moisture and chemicals. No drink which tests positively should be consumed. You must also be aware that The Drink Detective(TM) only tests for benzodiazepines, illicit GHB and ketamine. A drink, even one that has produced a negative test result, may still be unsafe as it may contain dangerous substances other than those tested for. Please take care. If you are in any doubt, do not drink. The manufacturer/seller does not make any representation or warranty in respect of this product save as specifically set out above. The manufacturer/seller cannot be held responsible for the presence of any drug other than as specified above. The manufacturer/seller excludes all liability for consequential and all other loss to the extent permitted by law. The manufacturer's/seller's liability for a defective product is limited to replacement or refund of the price paid.

Outside

Inside