

Document made available under the Patent Cooperation Treaty (PCT)

International application number: PCT/GB04/005261

International filing date: 16 December 2004 (16.12.2004)

Document type: Certified copy of priority document


Document details: Country/Office: GB
Number: 0405646.1
Filing date: 12 March 2004 (12.03.2004)

Date of receipt at the International Bureau: 24 January 2005 (24.01.2005)

Remark: Priority document submitted or transmitted to the International Bureau in compliance with Rule 17.1(a) or (b)



World Intellectual Property Organization (WIPO) - Geneva, Switzerland
Organisation Mondiale de la Propriété Intellectuelle (OMPI) - Genève, Suisse

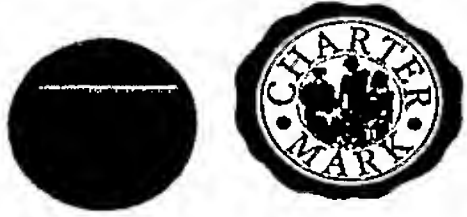

INVESTOR IN PEOPLE


The Patent Office
Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

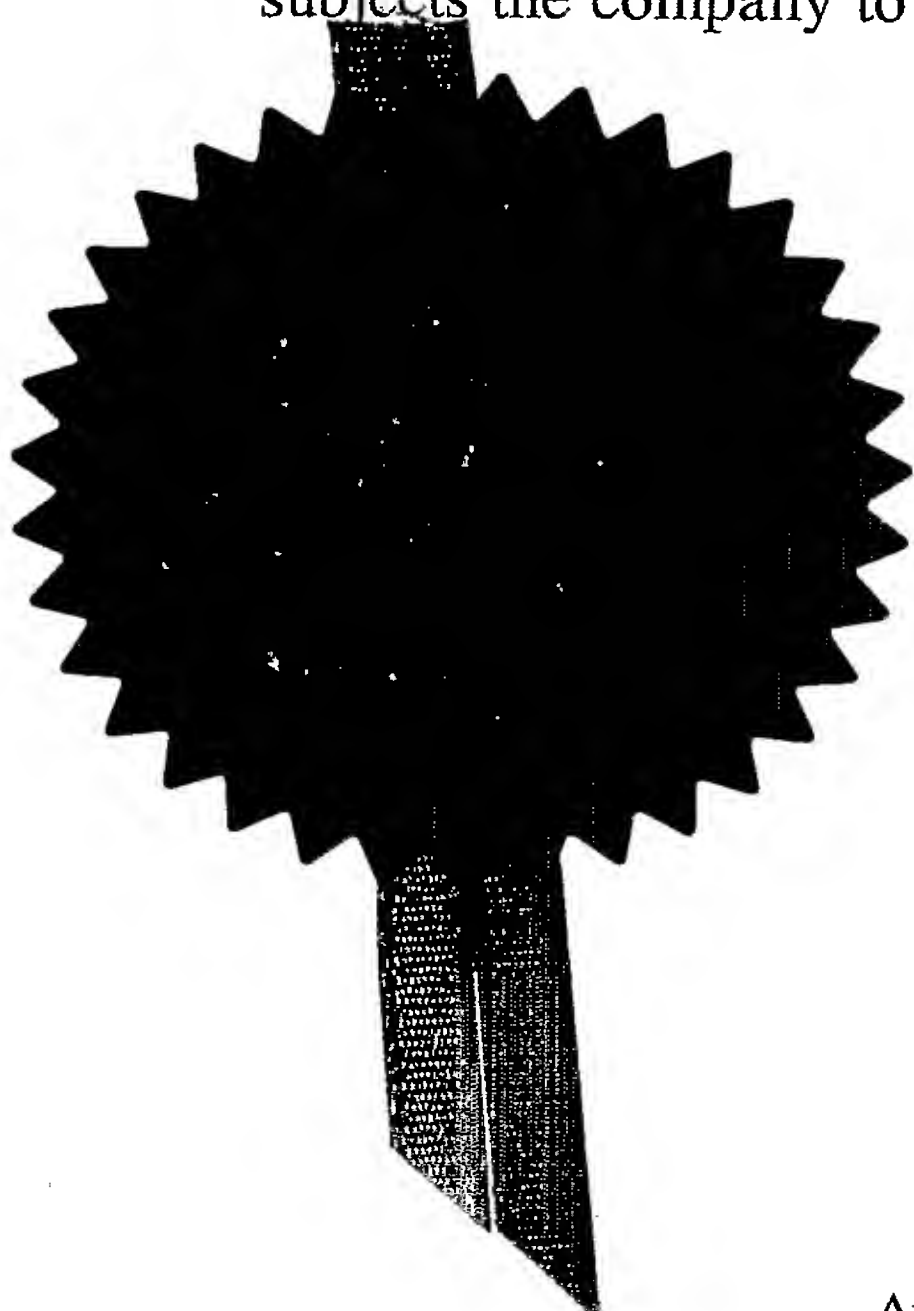

Signed

Dated 11 January 2005




Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

The Patent Office

Cardiff Road
Newport
South Wales NP10 8QQ

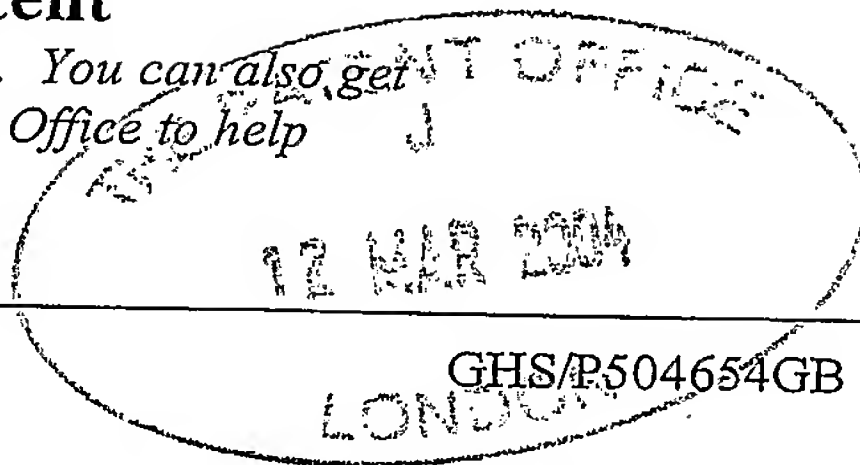

1.	Your reference	GHS/P504654GB		
2.	Patent application number (The Patent Office will fill in this part)	0405646.1		12 MAR 2004
3.	Full name, address and postcode of the or of each applicant (<u>underline all surnames</u>)	Bloomsbury Innovations Ltd 1 Prime Parkway Prime Enterprise Park Derby, Derbyshire DE1 3QB		
	Patents ADP number (if you know it)			
	If the applicant is a corporate body, give the country/state of its incorporation	Derbyshire		
4.	Title of the invention	Apparatus		
5.	Name of your agent (if you have one)	W.P.THOMPSON & CO.		
	"Address for service" in the United Kingdom to which all correspondence should be sent (including the postcode)	55 Drury Lane London WC2B 5SQ		
	Patents ADP number (if you know it)	158007		
6.	If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (if you know it) the or each application number	Country	Priority application number (if you know it)	Date of filing (Day/month/year)
7.	If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application	Number of earlier application		Date of filing (Day/month/year)
8.	Is a statement of inventorship and of right to grant of a patent required in support of this request? (Answer 'yes' if: a) any applicant named in part 3 is not an inventor, or b) there is an inventor who is not named as an applicant, or c) any named applicant is a corporate body. See note (d))	Yes		

8828362001

9. Enter then number of sheets for any of the following items you are filing with this form. Do not count copies of the same document.

Continuation sheets of this form	0
Description	9
Claims(s)	3
Abstract	0
Drawing(s)	4

10. If you are also filing any of the following, state how many against each item.

Priority documents

Translations of priority documents

Statement of inventorship and right to grant of a patent (Patents Form 7/77)

Request for preliminary examination and search (Patents Form 9/77)

Request for substantive examination (Patents Form 10/77)

Any other documents (Please specify)

11.

I/We request the grant of a patent on the basis of this application
 Signature W.P. Thompson & Co Date: March 12, 2004

W.P. THOMPSON & CO.

12.

Name and daytime telephone number of person to contact in the United Kingdom

Gill Smaggasgale
 020 7240 2220

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 01645 500505
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- For details of the fee and ways to pay please contact the Patent Office.

APPARATUS

The present invention relates to apparatus for use in testing the presence of foreign substances in a beverage. More particularly it relates to a test kit for detecting the presence of one or more
5 drugs in a beverage.

For ease of reference in the present application the term "drug" and "drugs" will be used to refer to any material which has been added to a beverage without the consent of the consumer. As will be explained in more detail below, the term will include alcohol, prescription pharmaceuticals,
10 over the counter pharmaceuticals, illegal substances and the like.

It has been known for many years for perpetrators wishing to cause harm to a victim to add a drug to the victim's beverage to alter the victim's behaviour or to incapacitate them totally. For example, alcohol may be added to a non-alcoholic drink or drugs may be added to either an
15 alcoholic or non-alcoholic drink. When the victim has drunk the beverage, their mood may be altered, or they may become totally incapacitated. While the victim is suffering from the effect of the drug the perpetrator may take advantage of the victim. In its mildest form, the victim's behaviour may simply be the source of amusement. Whilst this may be the perpetrator's intent, the action can have more serious consequences since the drug dose is generally unknown, as is
20 the effect of the combination of the drug with alcohol. Further the drug may have deleterious effects on any medical condition that the victim may have or may have an adverse interaction with any other medication that the victim may be taking. Unfortunately, there is usually a more sinister intent such as robbery, rape or even murder.

25 Drugs currently commonly used include alcohol, drugs from the benzodiazepine family, gamma hydroxybutyrate, 3,4-methylenedioxy-N-methylamphetamine (also known as MDMA or Ecstasy), ketamines, cannabis and the like. Recent studies suggest that there are at least two hundred drugs available legally or illegally which could be secretly applied to beverages for the purpose of disabling or incapacitating the victim.

30

Of particular concern at the present time is the use of flunitrazepam which is a potent benzodiazepine drug sold under the trade name Rohypnol. Whilst this drug has legitimate uses, it

has become known as the "date rape" drug because of its use as an additive to beverages to cause the victim to become willing to any suggestion and being unable to remember the situation afterwards. It is therefore commonly used to enable the perpetrator to obtain sex from the victim.

5 In an attempt to overcome this problem, the manufacturers of Rohypnol have added a blue dye to the prescription tablets to prevent their being added to drinks without the user's knowledge. However, there are still illicit forms of flunitrazepam available which are not coloured.

Drug rape usually goes unchecked since the victims are unable to remember the events or details
10 and therefore are unreliable witnesses in any legal action. Even if the victim can remember some details or have suspicions they are often unwilling to come forward.

It is therefore desirable to provide a testing means which will enable users to test their beverage to establish whether it is safe to drink. The benefits of this system include that the user knows
15 that he or she is safe. A further benefit of a testing method is that where the presence of a drug is detected, there is a possibility that the perpetrator may be identified. Any such device could also be used by authorities to screen drinks at the scene of an alleged offence.

Whilst there are a variety of established technologies available for testing for drugs, these are
20 generally based on laboratory or clinical tests and are therefore not suitable for use by the average consumer. They are also usually time consuming and expensive to run and require the use of clean rooms and for strict protocols to be observed. In addition, many of these systems are not able to function in the presence of an acid or alcohol and therefore have limited use in testing for drugs in beverages. These known tests also generally are only able to test for one drug at a time
25 and require resetting and recalibrating for the next drug to be tested.

Various proposals have been made to provide test apparatus which can be used at the site and time of consumption of the beverage. In US 6153147 a beverage analysis device is described in which a portion of liquid from the beverage can be captured and drawn into an analysis chamber
30 where it is subjected to a chemical reagent composition and the results of the colour assay are viewed through a window.

In US 2001/0046710 a test strip, toothpick or cotton swab is described which can be immersed into the beverage and which will exhibit a colour change when gamma hydroxybutyrate is present in the beverage.

- 5 In GB 0217339.1 filed 26th July 2002 a test strip is described which is particularly useful for the detection of gamma hydroxybutyrate or methylene dioxymethamphetamine. The test strip may be replaced with a coaster having a plurality of testing means.

10 Whilst these arrangements go some way to providing a suitable solution to the desirability of providing a test which can be used at the point of consumption, it is still desirable to provide alternative, and preferably improved, arrangements which preferably enable the user to test for one or more of a variety of drugs and which are preferably cost-effective to manufacture and easy to operate.

- 15 Thus according to the present invention there is provided a test kit for detecting the presence of one or more drugs in a beverage comprising a support having a plurality of detector strips attached thereto. The strips are preferably releasably attached thereto.

20 In a preferred arrangement of the present invention the kit will contain detector strips for detecting the presence of at least two, and preferably at least three, different drugs. Thus the kit may include at least one strip for detecting the presence of benzodiazepines and at least one detector strip for testing for gamma hydroxy butyrate. It may also include at least one strip for detecting the presence of alcohol.

- 25 In a further arrangement further strips may be included to test for one or more of ketamines, 3,4-methylenedioxy-N-methylamphetamine, cannabis, cocaine, barbiturates, opiates, tricyclic antidepressants, acetaminophen, propoxyphene and phencyclidine.

30 Whilst each detector strip may be configured to detect for more than one drug, in a preferred arrangement, each strip will test for one drug or one drug type only and will therefore be able to be formed of the optimum material for testing for the particular drug.

In one arrangement the kit will include at least two of each type of drug detector strips.

The detector strip may be of any suitable configuration. It will generally be of elongate configuration but other suitable configurations may be used. Where appropriate the strip may be rectangular or square. The strip will comprise the detection means. The strip may be completely supplied as the detection means or may comprise a bed supporting the detection means. The detection means will include assays or reagents which enable the target drug to be detected. For example, immunoassay, chemical spot, enzyme and the like techniques may be employed. For ease of reference, these will be collectively referred to as the "detection system".

The material from which the detection strip is constructed will depend on the detection system to be used. Where more than one detection system is used on different strips, they may be formed from different materials.

The antibodies used in immunoassay techniques are often expensive and have poor tolerance to acids, alcohol, moisture and temperature extremes. However, immunoassay tests may be developed which are suitable for use in the present invention.

In one arrangement of the present invention a detector strip may be present which is suitable for the detection of the presence of benzodiazepines. In one arrangement this may be prepared using an immunoassay technique for example a lateral flow competitive immunoassay or lateral flow direct binding immunoassay. Full details of one example of this immunoassay technique can be found in GB0217339 which is incorporated herein by reference.

In brief, a membrane is applied to the strip onto which is placed a pad of gold conjugate, bound to a first antibody having an affinity to the drug being tested, and a stripe, line or other indication from a second antibody is added which has an affinity for the first antibody. When the material comes into contact with the beverage the liquid will travel along the membrane carrying the gold conjugate-antibody complex with it as it migrates. If a drug is present, it will bind to this complex and block any binding sites on the complex. As the complex passes the location of the second antibody, the gold conjugate, which is magenta in colour, does not bind and passes beyond the line, so that no colour change is observed. The failure to produce a colour change would notify the user that there was a drug present.

If no drug is present, the gold conjugate-antibody complex does not have its binding sites blocked and as it passes the stripe of the second antibody binding occurs and a coloured stripe, line or other indication is obtained.

5

In another embodiment the production of binding to give a visual signal may be obtained by means of a sandwich assay instead of lateral flow technology. In this embodiment the constituents already described are placed through a porous membrane and the reaction takes place through the sandwich, revealing the result as a spot, line or other indication in the absence of a
10 drug in the benzodiazepines group. The technology may function either as a competitive sandwich immunoassay or as a direct sandwich immunoassay.

Examples of suitable "other indications" include a tick or a cross or a suitable word such as "pass" or "OK".

15

Immunoassay techniques may also be utilized to provide detector strips which may be used to detect for the presence of ketamines, 3,4-methylenedioxy-N-methylamphetamine, cocaine, barbiturates, opiates, tricyclic antidepressants, acetaminophen, propoxyphene and phencyclidine.

20 Whilst the detector strip for testing for benzodiazepines may be prepared using immunodiagnostic technology, other techniques for example colour change chemistry, may be used.

In the detector strips of the present invention, chemical spot techniques may be used as the
25 detection means. Chemical spot tests rely on a reaction, such as a redox reaction, between the substance being tested and the chemical used to perform the reaction which changes colour either because it produces a complex or because there is an alkalinity/acidity/pH change. Usually there is a colour change when the suspect drug is present. In one arrangement of the present invention the detector strip for detecting gamma hydroxybutyrate utilises chemical spot techniques.

30

In one embodiment of the present invention the test for gamma hydroxybutyrate operates by applying a small spot of chemical indicator on a bed. Any suitable material may be used for the

substrate including absorbent paper, cellulose sheet or film, cardboard or the like. The beverage is then brought into contact with the chemical indicator.

5 In one arrangement the detector strip for the detection of gamma hydroxybutyrate may be in the form of an elongate bed with a small area impregnated with the chemical indicator material. In an alternative embodiment the detector strip may be absorbent paper impregnated with the chemical indicator material.

10 In one arrangement of this latter embodiment the detector strip comprises a strip of absorbent paper impregnated with ferrous chloride in a buffer. When this strip is brought into contact with a beverage containing gamma hydroxybutyrate the ferrous chloride is oxidised to ferric chloride by the presence of the gamma hydroxybutyrate.

15 In another embodiment the chemical will react to pH change. This is possible since gamma hydroxybutyrate is generally a clandestine product derived from butyrolactone solvent with caustic, and is left in mildly alkaline form. Its addition to beverages will change the alkalinity of the beverage. This is readily detectable using the proposed chemical formulation. When there is typically 1 mg or more of gamma hydroxybutyrate in the beverage, the strip will turn blue when it comes into contact with the beverage.

20

This chemical spotting technology may also be used to provide a detection strip for use in the detection of amphetamines both pharmaceutical grade amphetamines and also those from 'street drugs' samples or clandestine tablets. In addition, it may be used to detect the presence of ketamines. Where a solution of Dragendorff reagent is impregnated in the absorbent material, it will produce a reddish colour when amphetamines are present or pink if ketamine is present. In addition, it will produce blackened speckles and tiny red/orange intense speckles when benzodiazepine is present. The blackened speckles are produced by the presence of bulking agents in the tablet and the red/orange intense speckles are due to the presence of the active benzodiazepine.

30

When aqueous solutions are applied to material impregnated with Dragendorff reagent blackening will occur although this blackening reduces where the solution is acidic. Thus in one

arrangement in which Dragendorff reagent is used in the detection strip, an acidic salt may be incorporated in the paper to minimise this effect when the strip is brought into contact with beverages.

- 5 Alcohol may be detected by an enzyme-based reaction in which chemical changes occur, leading to a colour change.

The support to which the detection strips are releasably attached may be formed of any suitable materials. Suitable materials include those which are lightweight and have sufficient inherent
10 strength not to become crushed in use and include cardboard and plastics material. The support material may be of any suitable configuration. Whilst it will generally be rectangular any shape may be used.

The detection strips may be attached to the support by any suitable means. In one arrangement at
15 least a portion of the detection strip may be formed integrally with the support and may be connected thereto by a frangible flange. Thus in use the detection strip may be removed from the apparatus by breaking the detection strip from the support.

In an alternative arrangement, the detection strips may be attached to the support by means of an
20 adhesive. Any suitable adhesive may be used. In use the desired strip can be removed from the support by means of breaking the adhesive bond. In an alternative arrangement the detection strip may be removed from the support by breaking the strip above the point where it is attached to the support by means of the adhesive. For example where the detection strip is made of impregnated paper or has a cardboard base it may simply be torn from the support. However,
25 where the detection strip has for example a plastics support, it may be provided with an area of weakness where breakage can occur. In an alternative arrangement, the adhesive may be a weak adhesive which will allow the detection strip to be removed from the support by, for example, peeling. Suitable adhesives include the repositional adhesive manufactured by 3M for use in Post-ItTM notes.

30

The test apparatus of the present invention may be of any suitable configuration. In one arrangement, the strips may extend from the support in a finger-like arrangement. However, in

an alternative arrangement, the strips will be laid across the face of the support. Where appropriate, the strips may be spaced in any suitable configuration on the support.

In a most preferred arrangement, the kit includes a cover for the support to protect the detection strips when they are not in use. The cover may be a wallet or box into which the support may be placed. In an alternative arrangement, the cover is integral with the support. Thus in one arrangement the support may be sized such that it may be folded into two parts such that when "closed" a first part may be laid over the strips which are attached to the second part. In use, the first part will be folded back to reveal the strips.

In another arrangement, there may be a third part which at least partially folds over the first part when it is in a closed position to help hold it in position.

Thus in one particularly preferred arrangement, the kit of the present invention will be configured to resemble a match book. The overall size of the kit when closed may be similar to a conventional match book or may be of the size of a business card.

Where the support is made from non-plastics material such as of cardboard, it may be coated on one or both surfaces with a water-resistant coating to provide protection.

The support may also be printed with instructions for use and/or advice and/or guidance relating to drugs. Advertising material may also be present.

Thus a particular benefit of the present invention is that the user has a discrete kit providing them with a wide range of testing options and which is simple to use and operate.

Various modifications of the apparatus may be provided. Since it is often dark in venues, the apparatus may include means to enable the user to find the appropriate strip, such as raised portions on the stick shaped to correspond to the drug to be tested. Additionally or alternatively, fluorescent markings may be used or a small light device may be included.

The present invention will now be described, by way of example, with reference to the

accompanying drawings:

Figure 1 illustrates one arrangement of the present invention;

5 Figure 2 illustrates a second arrangement of the present invention;

Figure 3 illustrates a still further arrangement of the present invention;

Figure 4 illustrates a product of the present invention; and

10

Figure 5 illustrates an alternative arrangement for a product of the present invention.

As illustrated in Figure 1, the first arrangement of the present invention comprises a support 1
15 and a plurality of tear-off detection strips 2.

Figure 2 illustrates a second arrangement where the detector strips 2 are attached by means of an adhesive to a face of the support 1.

20 Figure 3 illustrates a preferred arrangement where the detector strips 2 are located on the face of a support 1. The support has folds 3a and 3b which enables it to be folded over the strips to protect the strips and to provide a match-book like configuration. An example of the product is shown in Figure 4.

25 As illustrated in Figure 5 each detector strip may be accompanied by instructions for that test. In this embodiment the liquid is applied to the strip rather than the strip being dipped into the drink. The strips in this configuration are longitudinal for the test for benzodiazepines and rectangular for the tests for gammahydroxybutyrate and ketamines.

CLAIMS

1. A test kit for detecting the presence of one or more drugs in a beverage comprising a support having a plurality of detector strips attached thereto.
2. A test kit according to Claim 1 wherein the plurality of detector strips are releasably attached to the support.
3. A test kit according to Claim 1 or 2 wherein the kit contains detector strips for detecting the presence of at least two or at least three different drugs.
4. A test kit according to Claim 3 where the kit includes at least one strip for detecting the presence of benzodiazepines and at least one detector strip for testing for gamma hydroxy butyrate.
5. A test kit according to any one of Claims 3 or 4 wherein the kit includes at least one strip for detecting one or more of ketamines, 3,4-methylenedioxy-N-methylamphetamine, cocaine, barbiturates, opiates, tricyclic antidepressants, acetaminophen, propoxyphene and phencyclidine.
6. A test kit according to any one of Claims 3 to 5 additionally including at least one strip for detecting the presence of alcohol.
7. A test kit according to any one of Claims 1 to 6 wherein each detector strip is configured to detect for more than one drug.
8. A test kit according to any one of Claims 1 to 6 wherein each detector strip will is configured to test for one drug or one drug type only.
9. A test kit according to any one of Claims 1 to 8 wherein the kit includes at least two of each type of drug detector strips.

10. A test kit according to any one of Claims 1 to 9 wherein the detection means include assays or reagents which enable the target drug to be detected.
11. A test kit according to Claim 10 wherein the detection means include immunoassay, chemical spot, and/or enzyme techniques.
12. A test kit according to Claim 11 wherein the detector strip uses colour change chemistry.
13. A test kit according to any one of Claims 1 to 12 wherein the detector strip is supplied as the detection means.
14. A test kit according to any one of Claims 1 to 12 wherein the detector strip comprises a bed supporting the detection means.
15. A test kit according to any one of Claims 1 to 14 wherein the support to which the detection strips are releasably attached is formed from cardboard or plastics material.
16. A test kit according to any one of Claims 1 to 15 wherein the or each detection is formed integrally with the support and is connected thereto by a frangible flange
17. A test kit according to any one of Claims 1 to 16 wherein the or each detection strips is attached to the support by means of an adhesive.
18. A test kit according to any one of Claims 1 to 17 wherein the kit includes a cover for the support.
19. A test kit according to Claim 18 wherein the cover is a wallet or box into which the support is placed.
20. A test kit according to Claim 18 wherein the cover is integral with the support.
21. A test kit according to Claim 19 wherein the support is sized such that when it is folded

into two parts a first part is laid over the strips which are attached to the second part.

22. A test kit according to Claim 19 wherein a third part is present which at least partially folds over the first part when it is in a closed position to help hold the first part in position.
23. A test kit according to any one of Claims 1 to 22 wherein the kit is configured to resemble a match book.
24. A test kit according to any one of Claims 1 to 23 wherein the support is made from non-plastics material and is coated on one or both surfaces with a water-resistant coating to provide protection.
25. A test kit according to any one of Claims 1 to 24 wherein the support is printed with instructions for use and/or advice and/or guidance relating to drugs.

1/4

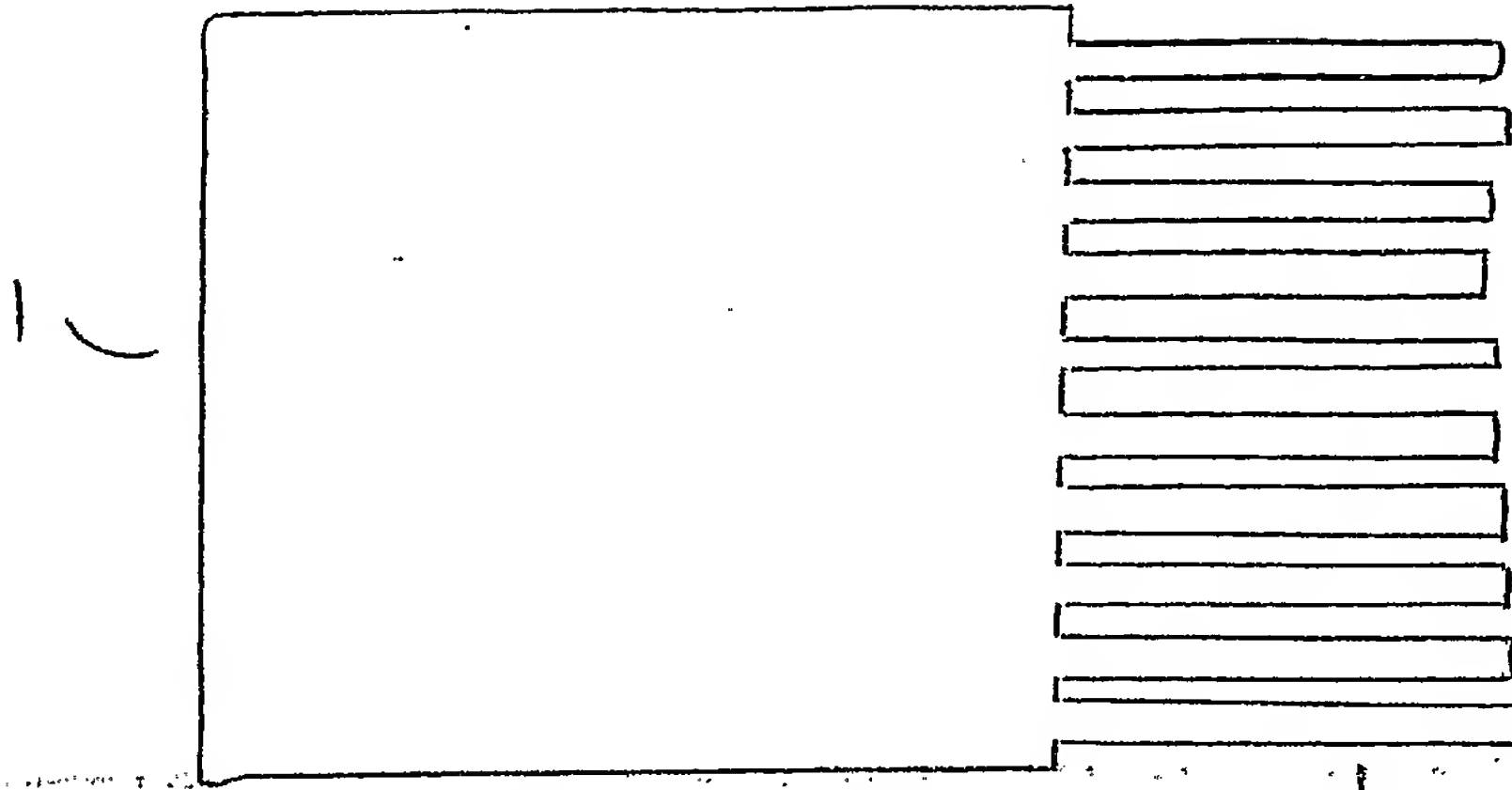


Figure 1

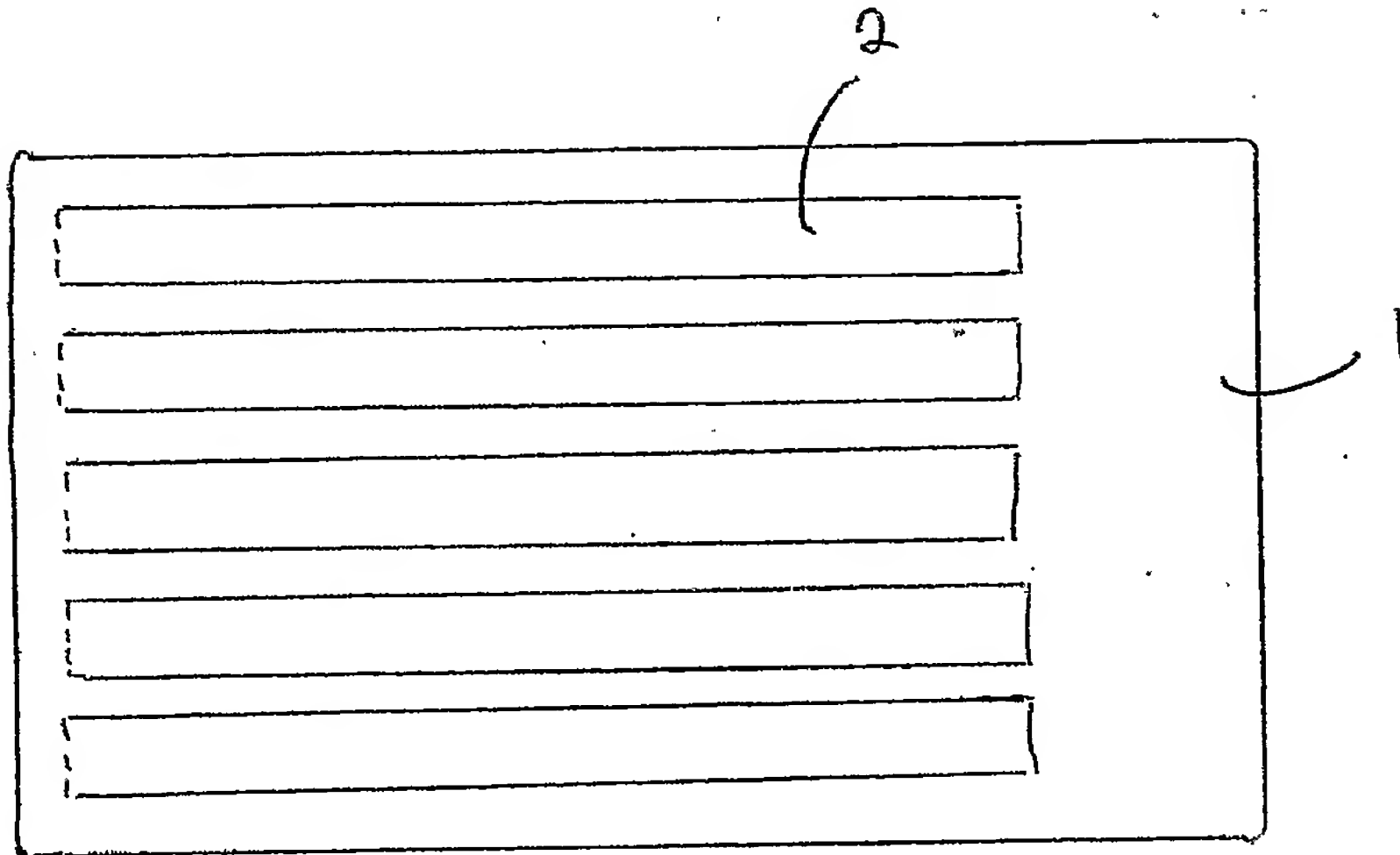


Figure 2




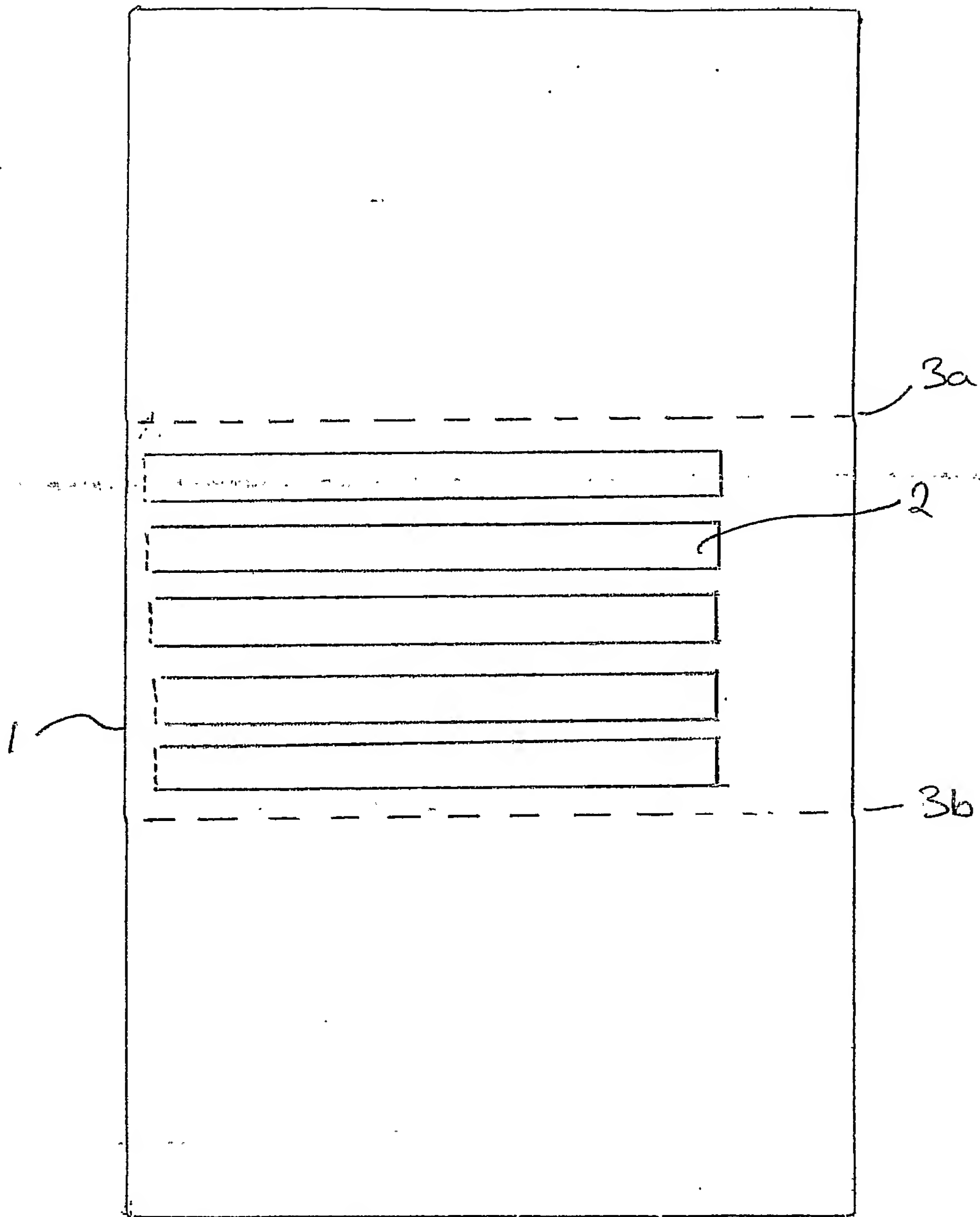



Figure 3





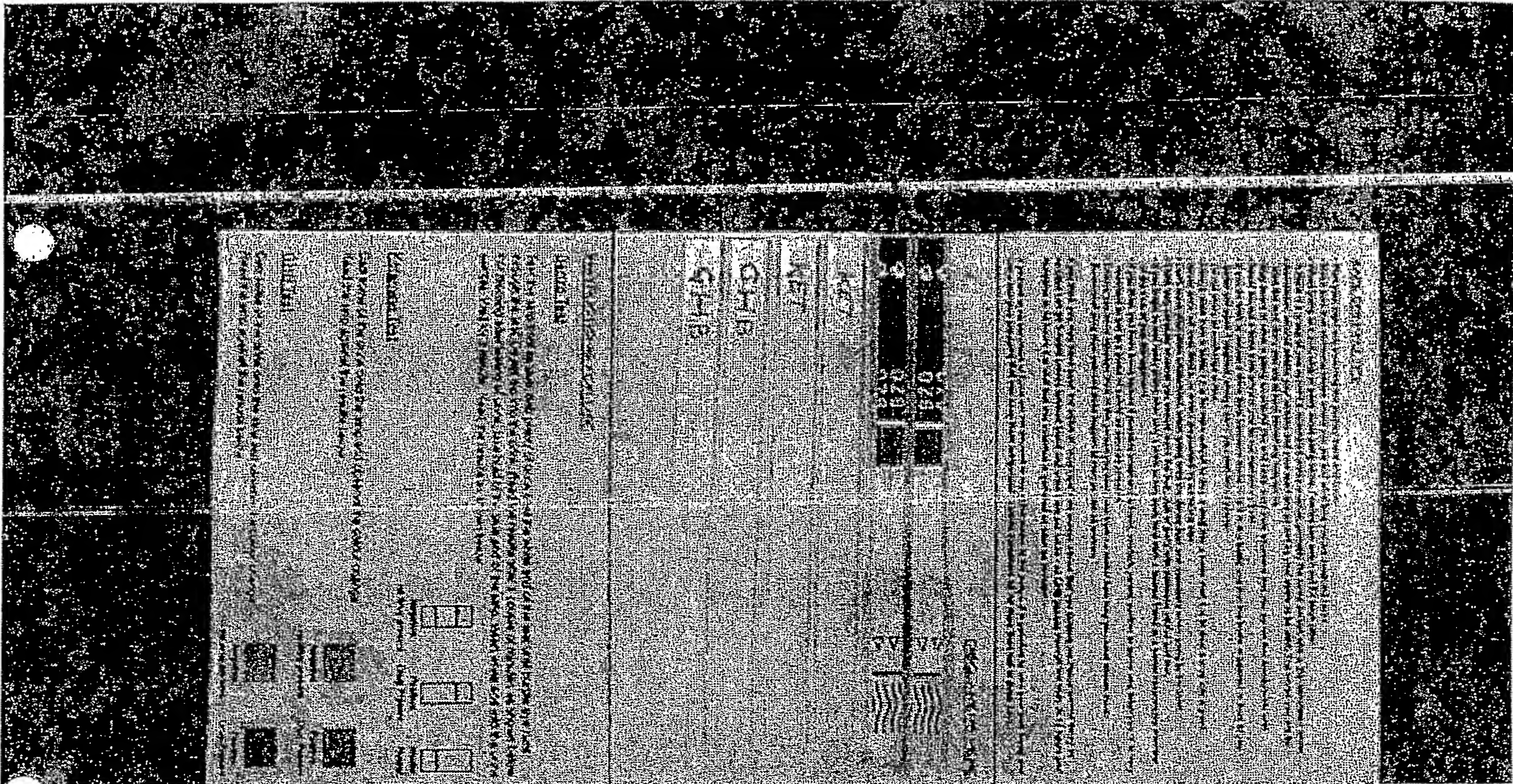

4/4

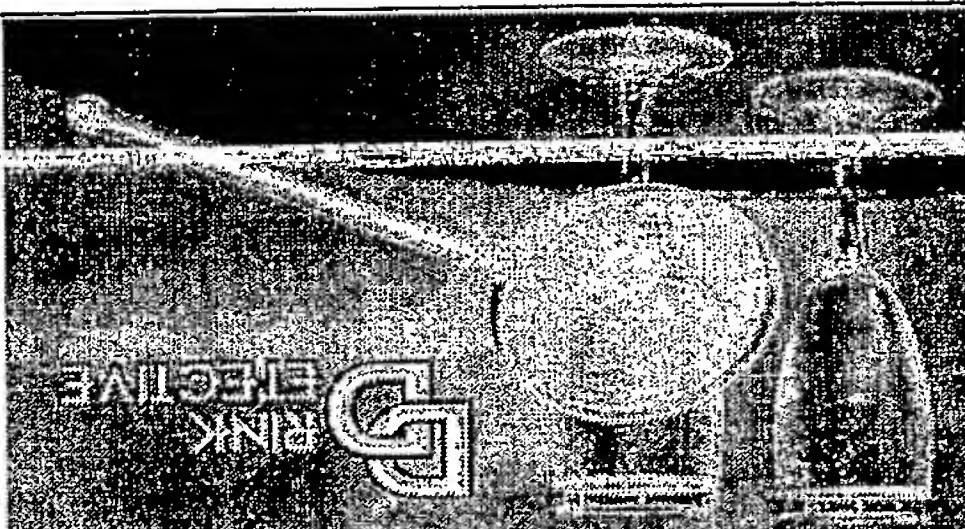

IF YOU HAVE A POSITIVE TEST RESULT

- DO NOT DRINK THE BEVERAGE. TRUST NO-ONE
- Give the drink to someone responsible. Consider phoning the police.
- A drug can take effect in as little as 20 minutes. If you drank from a spiked drink, ask a trusted friend to seek medical advice.
- If you are not with a trusted friend, tell the bar Manager and seek medical advice.
- IF YOU THINK YOU HAVE A SPIKED DRINK, YOUR ONLY CONCERN SHOULD BE YOUR PERSONAL SAFETY.

Sales help to support The Roofie Foundation's work with victims of drug-assisted assault. Visit www.roofie.com or call our helpline: +44(0)1332 362832

CAUTION
Instant Blue colour change

Put 1 drop of the drink onto the centre of the test pad.

HOW TO TEST

1 drop

Test here

Tear off strip goes here

RESULTS


© BILL 2004

OUTSIDE

HOW TO TEST

Put 1 drop of the drink onto the centre of the test pad.

K TEST

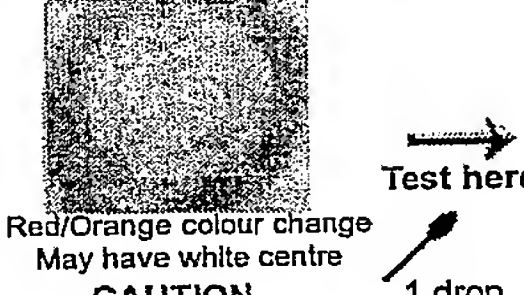

Red/Orange colour change
May have white centre

CAUTION

1 drop

Test here

Tear off strip goes here

HOW TO TEST

G TEST

Using the dropper, put 3 drops onto the test area as shown below (for best results leave 5 seconds between drops). After 2 minutes read against results shown below.

Tear off strip goes here

3 drops

Test here

RESULTS

2 lines - OK

1 line - CAUTION

WHAT THE DRINK DETECTIVE (TM) DOES AND ITS LIMITATIONS

The Drink Detective (tm) is designed to test for the presence of three substances: benzodiazepines (which include Rohypnol and Valium), illicit GHB (gamma-hydroxybutyrate) and ketamine. The Drink Detective (tm) does not test for any other drugs or substances. It is designed to test for minimal standard dosages of these substances and quantities below these levels may not be detected. Due to the sensitivity of the test, it is possible that the presence of certain substances could produce a false positive. In particular, it is possible that orange and tomato juices might produce a false positive for the 'K test'. The 'K test' might not pick up ketamine in drinks containing certain milk products. Milk products and soda waters containing bicarbonates of soda might produce a false positive for the 'G test'. Test your favourite drink and be alerted to different reactions. Using this test does not remove your responsibility to drink safely. Follow instructions exactly. Do not alter the test. Do not place in a drink or eat the test. Keep the test away from children, animals, moisture and chemicals. No drink which tests positively should be consumed. You must also be aware that The Drink Detective (tm) only tests for benzodiazepines, illicit GHB and ketamine. A drink, even one that has produced a negative test result, may still be unsafe as it may contain dangerous substances other than those tested for. Please take care. If you are in any doubt, do not drink. The manufacturer/seller does not make any representation or warranty in respect of this product save as specifically set out above. The manufacturer/seller cannot be held responsible for the presence of any drug other than as specified above. The manufacturer/seller excludes all liability for consequential and all other loss to the extent permitted by law. The manufacturer/seller's liability for a defective product is limited to replacement or refund of the price paid.

INSIDE

Figure 5

